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
1. Kerl JM, Ravenel JG, Nguyen SA, et al. Right heart: split-bolus injection of diluted contrast medium for visualization at coronary CT angiography. Radiology 2008;247:356-64.

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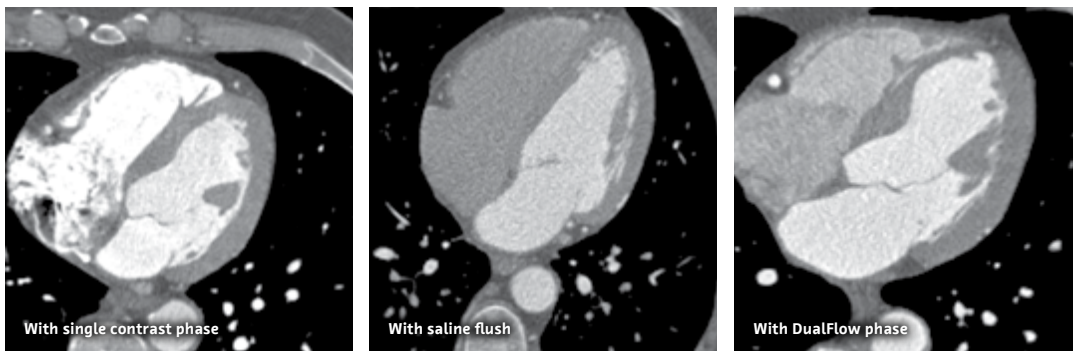Jensen, CT, "Dual Flow Contrast Injection for Coronary CTA Improves Visualization of the Right Heart,"

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
P3T® PA weight: 300 pound; contrast: 116 cc; main PA HU: 274

Standard weight 295 pound; contrast: 80 cc; main PA HU: 171

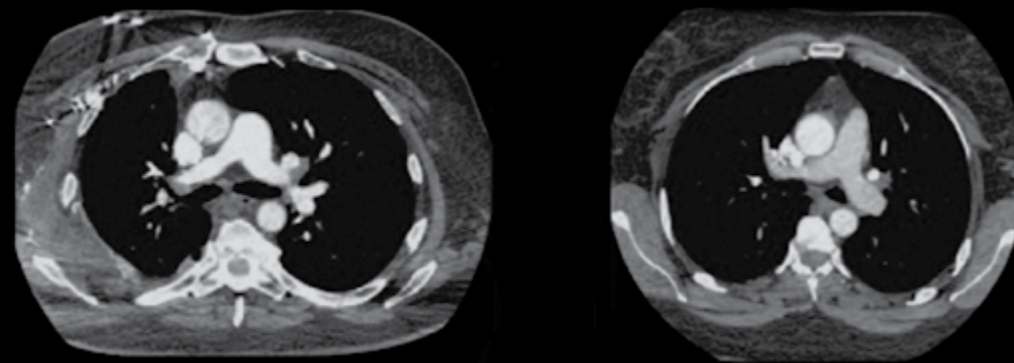When compared to standard protocols, P3T® PA protocols result in a higher percentage of exams ranked diagnostic without limitation.*

* Albeit at a slightly higher contrast dose than the standard protocol of 80 mL

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Item	Component	Material	Other Information
1	2 Multi-Patient Transfer Sets (MPTS)	Polyvinyl Chloride (PVC), Acrylonitrile butadiene styrene (ABS), Polystyrene (PS), Polycarbonate (PC), High Density Polyethylene (HDPE), Polypropylene (PP), Silicone Rubber	Use life: 12 hours Length: 115 cm Two check valves per MPTS
2	2 Syringe Assemblies (includes plunger)	Polyester (PET), PC, PP, Santoprene, Silicone	Use life: 12 hours Capacity: 200 mL
3	Expiration Label	Coated Vinyl	None
4	Multi-Patient Tubing Set (MPTS)	PVC, PC, PP	Use life: 12 hours
5	Single-Patient Tubing Set (SPTS)	PVC, Polycarbonate, Silicone Rubber, Silicone	Use life: 1 patient Length: 250 cm when fully extended Two check valves

General information	
Shelf Life: 4 years	Not made with Natural Rubber Latex
Sterilization Method: Radiation	SDS MP1 Primary Packaging: Polystyrene with Tyvek® cover
Pressure rated to 350 psi	SPD 250 Primary Packaging: Polyethylene with Tyvek® cover
Di(2-ethylhexyl)phthalate (DEHP) Free	Secondary Packaging: Cardboard Shipper
Non-Pyrogenic fluid path	

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Ordering Information

Systems listed include all standard cabling

Catalogue Number: SCT 211
MEDRAD® Stellant D Injection System
with Pedestal Head Mount

Catalogue Number: SCT 212
MEDRAD® Stellant D Injection System
with Counterpoise System

Catalogue Number: SCT 321*
MEDRAD® Stellant D CT Injection System
with Certegra® Workstation with a Pedestal Head Mount

Catalogue Number: SCT 322*
MEDRAD® Stellant D CT Injection System
with Certegra® Workstation with Counterpoise System Head Mount

Injection Specifications

Flow Rate (range & increments):
0.1–10 mL/sec in 0.1 mL Increments

Volume (range & increments):
1 mL to Syringe Capacity in 1 mL Increments

Programmable Pressure Limit (psi/kPa):
200 mL Syringe: 325 psi, 2241 kPa

Scan Delay:
0–300 Seconds (5 minutes) in 1 Second Increments

Pause:
1–900 Seconds (15 minutes) in 1 Second Increments

Hold:
Maximum HOLD Time is 20 Minutes

Syringes (Volume capacity):
200 mL Sterile Disposable Syringe

Maximum Number of Phases:
6

Miscellaneous Specifications

Electrical Requirements (VAC/Hz):
100–240 VAC, 50/60 Hz, 300 Volts-AMPs

Syringe Heat Maintainer:
35 degrees C +/-5 degrees, 95 degrees F +/-9 degrees

Head Mounting Options: Overhead Counterpoise System or Floor Pedestal (Floor Pedestal includes an Integral IV Pole)

DualFlow: Simultaneous Injection of Contrast and Saline is integrated in the CWS and included

Certegra® P3T® 2.0 Cardiac Application: Cardiac CTA Protocol Optimization Software

Certegra® P3T® 2.0 Abdomen Application: Abdominal Protocol Optimization Software

Certegra® P3T® 2.0 PA Application: Pulmonary Angiography Protocol Optimization Software

Informatics Starter Package: Access to Modality Worklist and Database of Injections

PACS Outbound Interface: Send Secondary Capture of Contrast Injection to PACS

RIS Outbound Interface: Send Contrast Injection Information to RIS

Speech Recognition Outbound Interface: Send Contrast Injection Information to Speech Recognition System

Imaging System Interfacing (ISI): Several Types of Interfacing with the Scanner

VirtualCare®: Remote service on the injector and software to ensure faster uptime